IMPORTANT SAFETY INSTRUCTIONS

DANGER – To reduce the risk of burns, electric shock and fire, this product must be used in accordance with the following instructions:

8. NEVER USE PAD WITHOUT REMOVABLE COVER IN PLACE.

9. O NOT USE IN AREAS WHERE OXYGEN IS BEING ADMINISTERED.

10. O NOT USE PINS OR OTHER METALLIC MEANS TO FASTEN THIS PAD IN PLACE.

11. DO NOT SIT ON OR CRUSH PAD – AVOID SHARP FOLDS.

12. EVER PULL THIS PAD BY THE SUPPLY CORD.

13. O NOT USE THE CORD AS A HANDLE.

14. AREFULLY EXAMINE INNER COVER AND LINE INNER COVERING OR LINE CORD SHOW ANY SIGN OF DETRIORATION.

15. DO NOT TAMPER WITH THIS PAD IN ANY WAY. THERE ARE NO USER-SERVICEABLE PARTS. IF SATISFACTORILY, RETURN IT TO AN AUTHORIZED SERVICE CENTER AT ADDRESS GIVEN.

16. O NOT USE THIS PAD WITH LINIMENT, SALVE OR OINTMENT PREPARATIONS THAT CONTAIN HEAT-PRODUCING INGREDIENTS, AS SKIN BURNS COULD RESULT. USE WHEREVER HOT APPLICATIONS ARE DESIRABLE FOR PERSONAL COMFORT AND WHEREVER RECOMMENDED BY YOUR PHYSICIAN OR DENTIST FOR THE RELIEF OF PAIN.

17. OOP CORD LOOSELY WHEN STORING; TIGHT WRAPPING MAY DAMAGE CORD AND INTERNAL PARTS.

18. K EP THE PAD AWAY FROM HEATED SURFACES.

19. LACE PAD ON TOP AND NOT UNDER THE PART OF THE BODY NEEDING HEAT.

20. O NOT PUT THIS PAD IN THE MICROWAVE.

21. O NOT LEAVE PAD UNATTENDED, ESPECIALLY IF CHILDREN ARE PRESENT.

22. O NOT USE PADS UNDERNEATH BEDCOVERS AND OTHER TYPES OF BLANKETS. DO NOT PLACE ABOVE PADS. ONLY THE CLOTH COVER PROVIDED WITH THE PAD SHOULD BE USED.

23. DO NOT USE PAD ON SENSITIVE SKIN OR WHEN TAKING MEDICATION THAT CAN CAUSE SENSITIVE SKIN.

24. AREFULLY EXAMINE INNER COVER AND LINE INNER COVERING OR LINE CORD SHOW ANY SIGN OF DETRIORATION.

25. DO NOT TAMPER WITH THIS PAD IN ANY WAY. THERE ARE NO USER-SERVICEABLE PARTS. IF SATISFACTORILY, RETURN IT TO AN AUTHORIZED SERVICE CENTER AT ADDRESS GIVEN.

26. O NOT USE THIS PAD WITH LINIMENT, SALVE OR OINTMENT PREPARATIONS THAT CONTAIN HEAT-PRODUCING INGREDIENTS, AS SKIN BURNS COULD RESULT. USE WHEREVER HOT APPLICATIONS ARE DESIRABLE FOR PERSONAL COMFORT AND WHEREVER RECOMMENDED BY YOUR PHYSICIAN OR DENTIST FOR THE RELIEF OF PAIN.

FEDERAL COMMUNICATIONS COMMISSION NOTICE

WARNING: Changes or modifications to this unit not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) This device must accept any interference received, including interference that may cause undesirable operation.

NOTE: The heat will automatically turn off after 60 minutes, at which point the green LED indicator light will flash continuously until the ON/OFF button is pressed.


NOTE: The heat will not exceed 30 minutes at one time. To turn off, press and hold ON/OFF button until all lights turn off. For maximum benefit, use should not exceed 30 minutes at one time.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications; however, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Standards: FCC Part 15 Subpart B, ANSI C63.4-2014
1-800-334-4031
Responsible Party: Conair Corporation, 1 Cummings Point Road, Stamford, CT 06902

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) This device must accept any interference received, including interference that may cause undesirable operation.

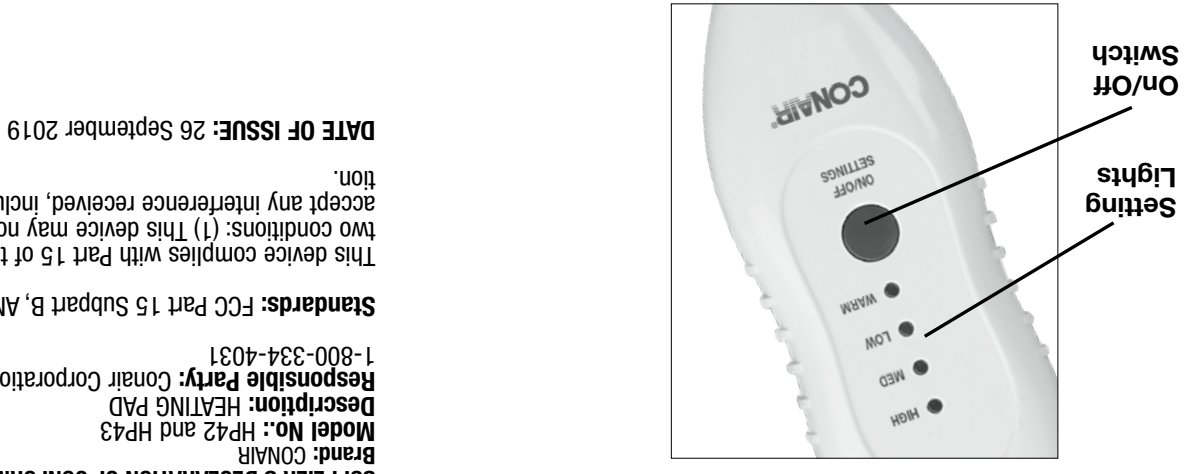
1. Plug the electrical cord into a 120V AC electrical outlet.
- NOTE: When first plugged into the main power supply, all lights will flash once. This to make sure the controller functions properly.
2. Push the ON/OFF button; the power indicator light will illuminate to signal that the heating pad is at Warm setting. Press the settings buttons to choose Warm, Low, Medium or High levels.
3. After choosing your heat setting, place heating pad on top of the area to be treated.
4. To turn off, press and hold ON/OFF button until all lights turn off. For maximum benefit, use should not exceed 30 minutes at one time.
- NOTE: The heat will not exceed 30 minutes at one time. To turn off, press and hold ON/OFF button until all lights turn off. For maximum benefit, use should not exceed 30 minutes at one time.
- NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications; however, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:
- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

MOIST HEAT TREATMENT

1. Ensure heating pad is unplugged.
2. Wet the absorbent inner cloth and squeeze out excess water.
3. Insert the damp pad between the heating pad and the cloth cover.
4. Plug the electrical cord into a 120V AC electrical outlet.
5. Push the ON/OFF button; the power indicator light will illuminate to signal that the heating pad is at Warm setting. Press the settings buttons to choose Warm, Low, Medium or High levels.
6. Place heating pad on top of the area to be treated.
7. When treatment is done, press and hold the ON/OFF button to turn unit off and unplug. For maximum benefit, use should not exceed 30 minutes at one time.
8. After use, dry the cloth cover, inner cloth and the heating pad before storing.

NOTE: This heating pad has a 60-minute auto-off feature. During use, check skin under heating pad frequently. Do not leave heating pad unattended.

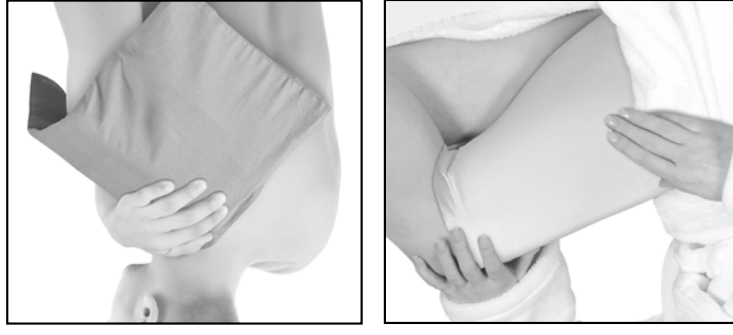
RESTRICTED HEATING PAD USES

For additional warnings, please see the Important Safety Instructions and Caution sections.

SAVE THESE INSTRUCTIONS

- Machine-ashable Cloth Cover (Figure 1)
- Heating Pad (Figure 2)
- Inner Cloth for Moist Heat (Figure 3)

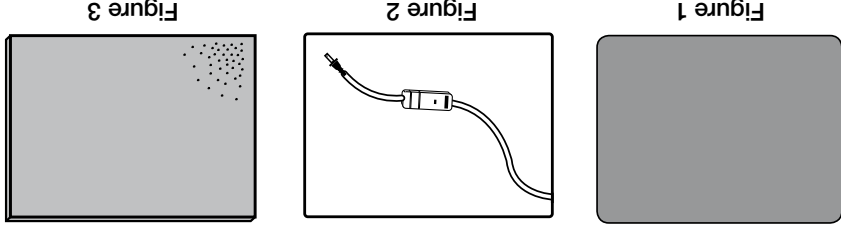Cautions:

- The cloth cover must always be used with the heating pad. An absorbent inner cloth is provided for moist heat application.
- This heating pad should not be used in the following cases: If you have diabetes, tuberculosis, benign or malignant tumors, phlebitis and thrombodosis, hemorrhages, inflammation of the skin, undetermined calf pain or joint pain, or if you have poor circulation. Consult a physician before using if you have any kind of fungal infection, swelling, fracture or persistent pain.

- Consult your physician should you have any therapeutic questions regarding use of this heating pad.
- Do not use on any unexplained pain, or swollen muscles, or following a serious muscular injury before consulting your physician.
- Keep away from children.
- Discontinue use of product if pain increases or swelling occurs after application.
- Do not use while bathing or showering. Use only in a dry area.
- Do not use your doctor before using in case of illness or medical condition, including but not limited to pacemaker use, pregnancy or cancer.
- If you experience pain in a muscle or joint for a prolonged period of time, discontinue use and consult your doctor. Persistent pain could be a symptom of a more serious condition.
- Never use any heating products on open wounds, discolored areas, or any area of the body that is swollen, burned, bruised, inflamed or where skin eruptions, ulcerated sores or varicose veins are present.
- This product should never be used by any individual suffering from any ailments that would limit or altogether inhibit the user's capacity to feel or have any sensation in the body part that this product is being used for.
- Use of this product should be pleasant and comfortable. Should pain or discomfort result, discontinue use and consult your physician.
- Ensure that your hands are dry when operating the controls or removing the plug.

LIMITED THREE-YEAR WARRANTY (U.S. AND CANADA ONLY)

Conair will repair or replace (at our option) your unit free of charge for 36 months from the date of purchase if the appliance is defective in workmanship or materials. To obtain service under this warranty, return the defective product to the service center listed, together with your purchase receipt. California residents need only provide proof of purchase and should call 1-800-3-CONAIR for shipping instructions. In the absence of a purchase receipt, the warranty period shall be 36 months from the date of manufacture.



ANY IMPLIED WARRANTIES, OBLIGATIONS, OR LIABILITIES, INCLUDING FITNESS FOR A PARTICULAR PURPOSE, SHALL BE LIMITED IN DURATION TO THE 36-MONTH DURATION OF THIS WRITTEN, LIMITED WARRANTY.

IN NO EVENT SHALL CONAIR BE LIABLE FOR ANY SPECIAL, INCIDENTAL, OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS OR IMPLIED. Some states do not allow the exclusion or limitation of special, incidental, or consequential damages, so the above limitation may not apply to you.

This warranty gives you specific legal rights, and you may also have other rights, which vary from state to state.

HOW TO ASSEMBLE

To assemble the heating pad, insert the heating pad completely into the cloth cover with the electric cord extending out from the cloth cover (see below). To secure in a desired position, attach each one end of the adjustable strap to the fabric cover. Wrap the fabric strap around desired location and attach the opposite end to the fabric cover and adjust as needed. If moist heat is desired, follow instructions in the Moist Heat Treatment section.

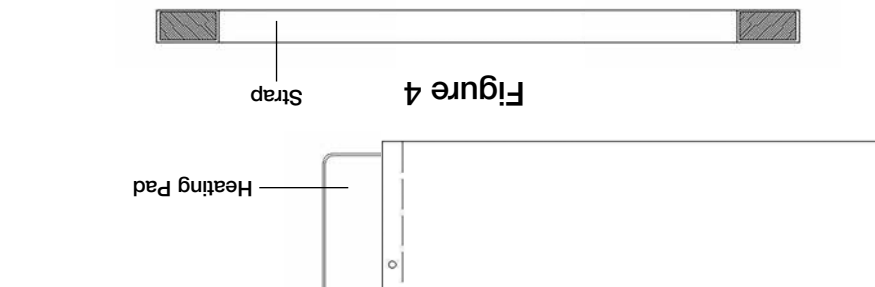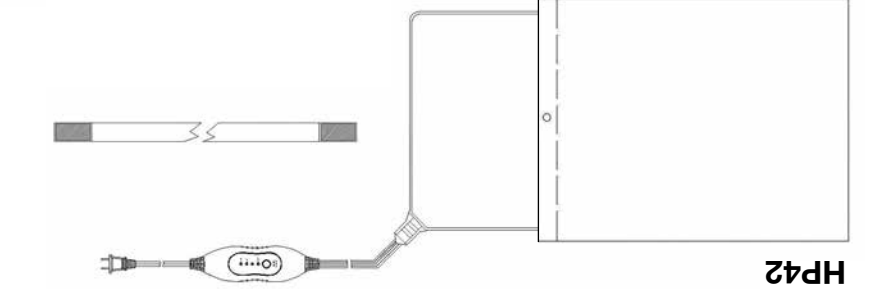
OPERATING INSTRUCTIONS

This heating pad comes assembled (Figure 4).

This appliance has a polarized plug (one blade is wider than the other). As a safety feature, this plug will fit in a polarized outlet only one way. If the plug does not fit fully in the outlet, reverse the plug. If it still does not fit, contact a qualified electrician. Do not attempt to defeat this safety feature. This heating pad is intended for household use only and it is designed to operate at 110-125 volts AC. Use on Alternating Current (60 Hertz) only.

CONAIR comfort™ deluxe/deluxe king size heating pad

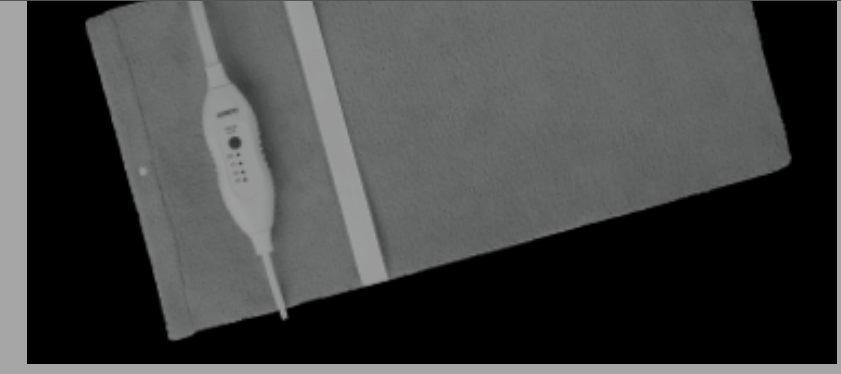Instruction booklet

For Models HP42 and HP43
Always read the instruction book carefully before using.

Service Center
7475 North Glen Harbor Blvd.
Glendale, AZ 85307
©2020 Conair Corporation
Please register this product at www.conair.com/registration
web at www.conair.com
For information on any Conair product, call 1-800-3-CONAIR or visit us on the
stock image
20PA073069
IB-16804

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