USER MAINTENANCE INSTRUCTIONS

CARING FOR YOUR TRUE GLOW® MICRODERMABRASION UNIT

1. After each use, wipe the device with a damp cloth to remove any oils or lotions used during treatment. Be sure you keep the charging contacts on the Microdermabrasion Unit clean.

2. After each use, clean the device to prevent buildup of oils and lotions. The 2 pin contacts must be clean for the unit to charge properly.

3. Clean the attachment head and unit after each use for sanitation purposes by rinsing under warm, soapy water to help remove oil and lotion buildup. Clean the attachment head and unit after each use with rubbing alcohol or antibacterial soap. Wipe with a dry cloth and store unit in its box or in a clean, dry place.

4. Do not put unit in dishwasher for cleaning.

TROUBLE-SHOOTING GUIDE

<p>THE MICRODERMABRASION UNIT IS NOT CHARGING.</p> <p>NO WHITE FLASHING CHARGING LIGHT.</p> <p>THE CHARGING LIGHT DOES NOT CHANGE TO SOLID WHITE.</p>	<ol style="list-style-type: none"> 1. Make sure the contacts on the unit and charging cable are clean. See User Maintenance Instructions for cleaning procedure. 2. Make sure the unit is plugged in properly. 3. Check the outlet and make sure the power is on. Most bathroom outlets have a test and reset button. Press the reset button and recheck the outlet.
<p>THE MICRODERMABRASION UNIT WILL NOT TURN OFF.</p>	<p>You must hold the ON/OFF button for more than 1 second to turn the unit off.</p>
<p>THE MICRODERMABRASION UNIT OR ATTACHMENTS FALL OFF DURING USE.</p>	<p>Make sure the attachments are snapped on properly.</p>

IMPORTANT SAFETY INSTRUCTIONS

When using electrical appliances, especially when children are present, basic safety precautions should always be taken, including the following:

READ ALL INSTRUCTIONS BEFORE USING THE MICRODERMABRASION UNIT

DANGER — Any appliance is electrically live, even when the switch is off.

WARNING — To reduce the risk of burns, electrocution, fire, injury or death to persons:

1. Always unplug appliance from electrical outlet immediately after using (except when charging unit). Do not allow to be used as a toy.
2. Do not place or store appliance where it could fall or be pulled into a tub or sink.
3. Do not use appliance while bathing or in the shower.
4. Do not place appliance in, or drop into water or other liquids.
5. If appliance has fallen into water, do not reach for it. Unplug immediately.
6. Unplug this appliance before cleaning.
7. Do not operate the charger if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged or dropped into water.
8. Do not place or store the charging cable where it can fall or be pulled into a tub or sink.
9. Do not immerse charging cable in water or any other liquids.
10. If charging cable falls into water or another liquid, unplug immediately.
12. Not for use by children.

13. Do not pull, twist or wrap the cord around the appliance.
14. Never drop or insert any object into any opening on the appliance.
15. THIS PRODUCT IS FOR HOUSEHOLD USE ONLY.
16. Do not leave appliance when plugged in. Unplug from outlet when not in use and before cleaning and servicing.
17. Do not unplug by pulling on cord. To unplug, grasp the plug, not the cord.
18. Do not handle plug or appliance with wet hands.

CAUTION —

1. Dry location use only. Do not expose to liquids, vapor, or rain.
2. Ensure jack is plugged into the unit before plugging the USB cable into the outlet.
3. Do not operate any appliance that has a damaged output cord or housing.
4. Always unplug the unit before cleaning or when not in use.
5. Be sure voltage is correct before plugging in. The rated input is 5V DC 200mA.

The True Glow® Microdermabrasion Unit handle, charger and cable have no serviceable parts. If the charger cord is damaged, it must be unplugged and replaced. Call 1-800-3-CONAIR for instructions or to purchase replacement brush heads.

SAVE THESE INSTRUCTIONS

SENSITIVITY AND CAUTIONS SPECIAL CAUTIONS FOR DELICATE AREAS

The True Glow® Microdermabrasion Unit is intended to gently remove dirt and makeup from your pores without irritation, resulting in healthy-looking skin.

BATTERY REMOVAL

1. Using a small Phillips-head screwdriver, remove the screws at the front of the unit.
2. Remove the control panel and cut the battery out. Bring the battery to the nearest collection center for proper disposal. DO NOT throw the battery into the general municipal waste system.

CONTAINS LITHIUM ION BATTERY. MUST BE COLLECTED, RECYCLED OR DISPOSED OF PROPERLY. DO NOT DISPOSE OF IN THE MUNICIPAL WASTE. KEEP AWAY FROM WATER

ALWAYS UNPLUG SMALL APPLIANCES

RBRC™ is a trademark of the Rechargeable Battery Recycling Corporation.

The batteries are sealed inside the True Glow® Microdermabrasion Unit, so they cannot be easily removed for recycling. The entire True Glow® Microdermabrasion Unit should be returned when recycling the batteries for this product.

LIMITED ONE YEAR WARRANTY (U.S. AND CANADA ONLY)

Conair will repair or replace (at our option) your unit free of charge for 12 months from the date of purchase if the appliance is defective in workmanship or materials.

To obtain service under this warranty, return the defective product to the service center listed on back that is nearest you, together with your purchase receipt. California residents need only provide proof of purchase and should call 1-800-3-CONAIR for shipping instructions. In the absence of a purchase receipt, the warranty period shall be 12 months from the date of manufacture.

ANY IMPLIED WARRANTIES, OBLIGATIONS, OR LIABILITIES, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, SHALL BE LIMITED IN DURATION TO THE 12-MONTH DURATION OF THIS WRITTEN, LIMITED WARRANTY.

Some states do not allow limitations on how long an implied warranty lasts, so the above limitations may not apply to you.

IN NO EVENT SHALL CONAIR BE LIABLE FOR ANY SPECIAL, INCIDENTAL, OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS OR IMPLIED, WHATSOEVER.

Some states do not allow the exclusion or limitation of special, incidental, or consequential damages, so the above limitation may not apply to you.

This warranty gives you specific legal rights, and you may also have other rights, which vary from state to state.

FEATURES

THREE SPEED OPTIONS

The three speeds available on your True Glow® Microdermabrasion Unit allow you to customize your skincare routine to suit your individual needs.

ATTACHMENTS

FINE ATTACHMENT: For sensitive skin.
COARSE ATTACHMENT: For deeper exfoliation.
PORE ATTACHMENT: For blackhead extraction.

NOTE: Your True Glow® Microdermabrasion Unit may only include certain attachments. See retailer for replacement units or call 1-800-3-CONAIR or visit www.trueglowbyconair.com

RECHARGEABLE

The True Glow® Microdermabrasion Unit is cordless and requires a full charge before use. After unit is fully charged, it can provide up to 35 minutes of use without recharging.


OPERATING PRECAUTIONS

There is a built-in, rechargeable Li-polymer battery. Do not throw into fire, apply heat, use or leave in a high-temperature environment.

CAUTION: Any changes or modifications not expressly approved by the Conair Corporation, the manufacturer of this product and the party responsible for compliance, could void the user's authority to operate the equipment.

WARNING: IF YOU THINK THE POWER IS OFF WHEN THE SWITCH IS OFF, YOU'RE WRONG.

KEEP AWAY FROM WATER

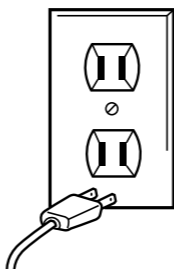Everyone knows that electricity and water are a dangerous combination. But did you know that an electrical appliance is still electrically live even if the switch is off? If the plug is in, the power is on. So when you are not using your appliances, keep them unplugged.

ALWAYS UNPLUG SMALL APPLIANCES

CHARGING THE MICRODERMABRASION UNIT

1. Before using the Microdermabrasion Unit for the first time, it must be fully charged for 2.5–3.5 hours. Plug the power cable into a wall outlet and plug the DC jack into the back of the unit. The white LED lights up, indicating that the unit is charging.

2. The charging indicator under the control panel will flash BRIGHT WHITE when charging. Once fully charged, the light will illuminate SOLID WHITE. Always check to ensure that the charging light illuminates when you plug the charging cable into the Microdermabrasion unit..

3. The WHITE light will start flashing for 5 seconds and then shut off when the unit reaches low battery capacity during operation.

4. The adapter may become warm while charging; this is normal operation.

5. In order to maintain the optimum capacity of the rechargeable battery, the Microdermabrasion Unit should be given a full, 2.5–3.5-hour charge every 3 months.

6. Unplug the adapter from the wall and the cable when not in use. Do not unplug the adapter from the wall when it is still connected to the unit. This will cause the battery to discharge

FCC MARKINGS

WARNING: Changes or modifications to this unit not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) This device must accept any interference received, including interference that may cause undesired operation.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

CAUTION: Any changes or modifications not expressly approved by Conair Corporation, the manufacturer of this product and the party responsible for compliance, could void the user's authority to operate the equipment.

STARTING THE MICRODERMABRASION UNIT AND SELECTING THE SPEEDS

Press the ON/OFF button 1 second and it will start at LOW setting, 2 times for MEDIUM setting and 3 times for HIGH setting. The white indicator will illuminate with different brightness (Low: dark; Med: bright; High: more bright) indicating the setting on the control panel.

STOPPING THE MICRODERMABRASION UNIT

At any time, the Microdermabrasion Unit can be stopped by pressing and holding the ON/OFF button for more than 1 second.

ATTACHING AND REMOVING MICRODERMABRASION HEADS

To attach the different attachments, hold the unit and push in the attachment onto the post until it snaps in place. To remove attachment heads, simply pull attachment off of the unit.

It is recommended to clean attachment head and unit after each use. If not clean, performance could be affected. Change the filter if the suction strength is weak.

SUPPLIER'S DECLARATION OF CONFORMITY

Brand: TRUE GLOW

Model No.: MD2

Description: Microdermabrasion Unit
Responsible Party: Conair Corporation, 1 Cummings Point Road, Stamford, CT 06902 (203) 351-9000

Standards: FCC Part 15 Subpart B, ANSI C63.4-2014

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

(1) This device may not cause harmful interference, and (2) This device must accept any interference received, including interference that may cause undesirable operation.

DATE OF ISSUE: TBD

Please register this product at www.conair.com/registration

SERVICE CENTER CONAIR CORPORATION

Service Department
7475 North Glen Harbor Blvd.
Glendale, AZ 85307

Visit us on the web at www.conair.com


©2020 Conair Corporation
Glendale, AZ 85307

stock.image

MD2 Models
20PS073376

IB-14717A

customize your treatment at home

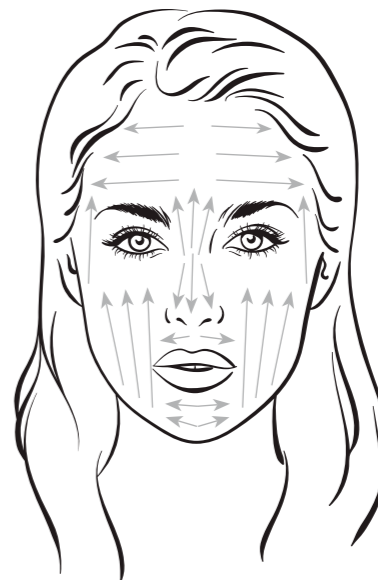OPERATING INSTRUCTIONS

1. Thoroughly cleanse and dry the skin.
2. Test on a small area of the neck or forearm.
3. Pull the skin taut with your free hand to allow the device to smoothly glide over the skin.
4. Move device in an upward, fluid motion across the skin. Do not hover in any spots.
5. Do not use directly above or below eyes.
6. Immediately after treatment, wash and dry the skin.
7. Apply a calming toner and recovery moisturizer.
8. Clean the attachment head and unit after each use with rubbing alcohol or antibacterial soap.
8. Wait 6 to 7 days between treatments.

NOTE: Do not do more than two passes over the same area during treatment.

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