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CONSERVAR ESTAS INSTRUCCIONES

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HOW TO ASSEMBLE

To assemble the heating pad, insert the heating pad completely into the cloth cover with the electric cord extending out from the cloth cover (Figure 4). If moist heat is desired, follow instructions in the Moist Heat Treatment section.

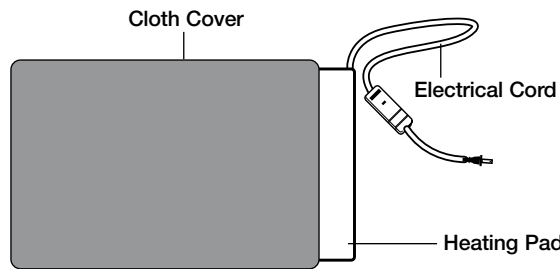
Figure 4

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Cautions:

- The cloth cover must always be used with the heating pad. An absorbent cloth is provided for moist heat application.
- This heating pad should not be used in the following cases: If you have diabetes, tuberculosis, benign or malignant tumors, phlebitis and thrombosis, hemorrhages, inflammation of the skin, undetermined calf pain or joint pain, or if you have poor circulation. Consult a physician before using if you have any kind of fungal infection, swelling, fracture or persistent pain.
- **Consult your physician should you have any therapeutic questions regarding use of this heating pad.**
- Do not use on any unexplained pain, or swollen muscles, or following a serious muscular injury before consulting your physician.
- **Keep away from children.**
- Discontinue use of product if pain increases or swelling occurs after application.
- Do not use while bathing or showering. Use only in a dry area.
- Consult your doctor before using in case of illness or medical condition, including but not limited to pacemaker use, pregnancy or cancer.
- If you experience pain in a muscle or joint for a prolonged period of time, discontinue use and consult your doctor. Persistent pain could be a symptom of a more serious condition.
- Never use any heating products on open wounds, discolored areas, or any area of the body that is swollen, burned, bruised, inflamed or where skin eruptions, ulcerated sores or varicose veins are present.
- This product should never be used by any individual suffering from any ailments that would limit or altogether inhibit the user’s capacity to feel or have any sensation in the body part that this product is being used for.
- Use of this product should be pleasant and comfortable. Should pain or discomfort result, discontinue use and consult your physician.
- Ensure that your hands are dry when operating the controls or removing the plug.

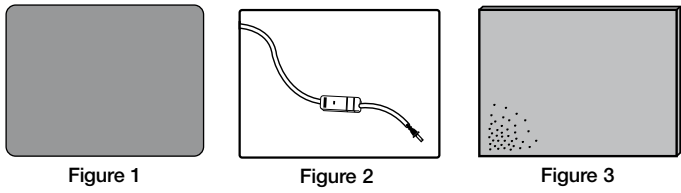
Figure 1

Figure 2

Figure 3

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

IMPORTANT SAFETY INSTRUCTIONS

DO NOT DESTROY

DANGER – To reduce the risk of burns, electric shock and fire, this product must be used in accordance with the following instructions:

1. Read all instructions.
2. **ALWAYS UNPLUG THE UNIT FROM THE ELECTRICAL OUTLET IMMEDIATELY AFTER USING.**
3. Do not use while sleeping. Avoid falling asleep while using the product.
4. Do not use on an infant or an animal.
5. This pad is not to be used on or by any individual with certain disabilities, sleeping or unconscious person, or a person with poor blood circulation unless carefully attended.
6. Do not use on areas of insensitive skin.
7. Burns can occur regardless of control setting; check skin under pad frequently.
8. Never use pad without removable cover in place.
9. Do not use in areas where oxygen is being administered.
10. Do not use pins or other metallic means to fasten this pad in place.
11. Do not sit on or crush pad – avoid sharp folds.
12. Never pull this pad by the supply cord.
13. Do not use the cord as a handle.
14. Carefully examine inner cover and line cord before each use. Discard the pad if inner covering or line cord show any sign of deterioration.
15. Do not tamper with this pad in any way. There are no user serviceable parts. If, for any reason, this pad does not function satisfactorily, return to authorized service center at address given.

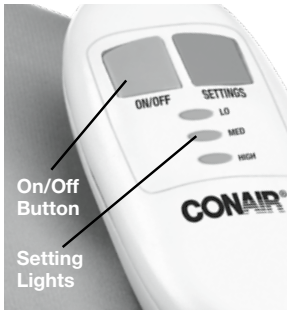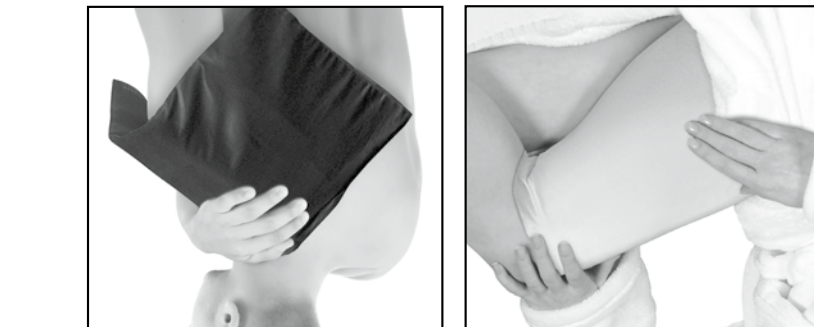Note: This heating pad has a 60-minute auto-off feature. During use, check skin under heating pad frequently. Do not leave heating pad unattended. After use, dry the cloth cover, absorbent cloth and the heating pad before storing.

8. After use, dry the cloth cover, absorbent cloth and the heating pad before storing.
Note: This unit has a 60-minute auto-off feature. For maximum benefit, use should not exceed 30 minutes at one time.

1. Ensure heating pad is unplugged.
2. Wet the absorbent cloth and squeeze out excess water.
3. Insert the damp cloth between the heating pad and the cloth cover.
4. Plug the electrical cord into a 120V AC electrical outlet.
5. Push the ON/OFF button, the power indicator light will illuminate to signal that the heating pad is at Low setting. Press the settings buttons to choose among Low, Medium and High levels.
6. Place heating pad on top of the area to be treated.
7. When treatment is done, press the ON/OFF button to turn unit off and unplug.

MOST HEAT TREATMENT

RESTRICTED HEATING PAD USES

Do not use while sleeping
 Do not use heating pad on children
 Do not use heating pad on chest

Do not sit on heating pad
 Do not sharply fold or crease heating pad
 Do not use under covers

For additional warnings, please see the Important Safety Instructions and Caution sections.

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FEDERAL COMMUNICATIONS COMMISSION NOTICE

WARNING: Changes or modifications to this unit not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment. This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) This device must accept any interference received, including interference that may cause undesired operation.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

CAUTION: Any changes or modifications to this unit not expressly approved by Conair LLC, the manufacturer of this product and the party responsible for compliance, could void the user's authority to operate the equipment.

SUPPLIER'S DECLARATION OF CONFORMITY

Brand: CONAIR®
Model No.: HP15X, ALL MODELS
Description: HEATING PAD
Responsible Party: Conair LLC, 1 Cummings Point Rd. Stamford, CT 06902 (203) 351-9000
Standards: FCC Part 15 Subpart B, ANSI C63.4-2014
 This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) This device must accept any interference received, including interference that may cause undesirable operation.
DATE OF ISSUE: XX/XX/XXXX

LIMITED THREE YEAR WARRANTY (U.S. AND CANADA ONLY)


Conair will repair or replace (at our option) your unit free of charge for 36 months from the date of purchase if the appliance is defective in workmanship or materials.

ANY IMPLIED WARRANTIES, OBLIGATIONS, OR LIABILITIES, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, SHALL BE LIMITED IN DURATION TO THE 36 MONTH DURATION OF THIS WRITTEN, LIMITED WARRANTY. Some states do not allow limitations on how long an implied warranty lasts, so the above limitations may not apply to you.

IN NO EVENT SHALL CONAIR BE LIABLE FOR ANY SPECIAL, INCIDENTAL, OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS OR IMPLIED, WHATSOEVER. Some states do not allow the exclusion or limitation of special, incidental, or consequential damages, so the above limitation may not apply to you.

This warranty gives you specific legal rights, and you may also have other rights, which vary from state to state.
 Conair LLC
 Service Center
 7475 North Glen Harbor Blvd.
 Glendale, AZ 85307
 Visit us on the web at www.conair.com
 ©2022 Conair LLC
 Glendale, AZ 85307
 stock.image
 21PA082319
 IB-17557

CONAIRcomfort™ heating pad

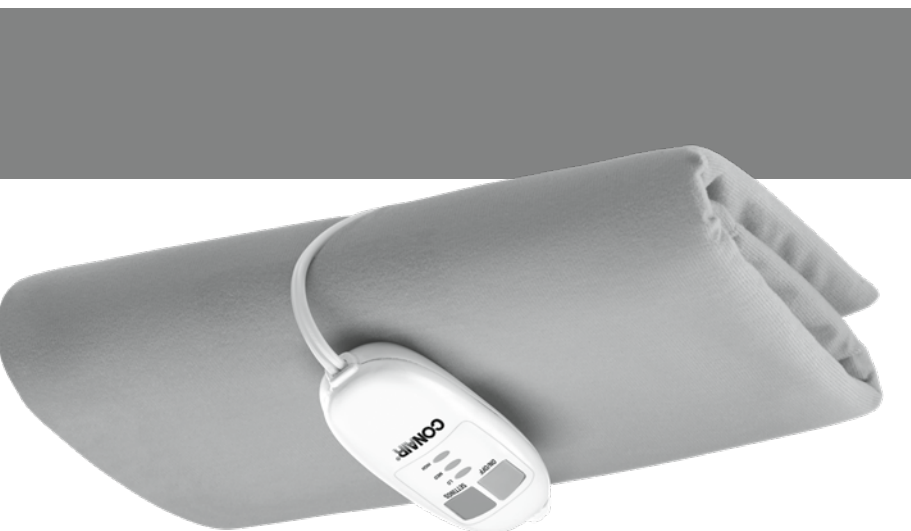for all HP15X models
 always read the instruction book carefully before using.
 For your safety and continued enjoyment of this product,
 21PA082319
 stock.image

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Centro de servicio
Conair LLC
Service Department
7475 North Glen Harbor Blvd.
Glendale, AZ 85307

Visítenos en www.conair.com

©2022 Conair LLC
 Glendale, AZ 85307

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IB-17557