Instruction Manual • Model WJ12 (All Versions)

For your safety and continued enjoyment of this product, always read the instruction book carefully before using.

Service Center:
Conair Corporation, Service Department, 7475 N. Glen Harbor Blvd., Glendale, AZ 85307
IB-1661
19105634
©2019 Conair Corporation, 7475 N. Glen Harbor Blvd., Glendale, AZ 85307
For more information regarding this product or call 1-800-334-4031
visit www.interplak.com

SAVE THESE INSTRUCTIONS

SAVE THESE INSTRUCTIONS

INTERPLAK[®]
cordless
PORTABLE
water flossing system

Limited 1-Year Warranty (U.S. and Canada Only)
Conair will repair or replace at our option, your unit free of charge for 12 months
from the date of purchase if the unit is defective in workmanship or materials.


MEDICAL WARNING
For your safety:

IMPORTANT SAFEGUARDS BEFORE USING READ ALL INSTRUCTIONS

- DANGER** • DO NOT submerge the unit in water or any other liquids. Use this product only for its intended use as described in this manual. Do not use it around children unless recommended by the manufacturer.
- WARNING** • Close supervision is necessary when this product is used by, on, or near children or individuals with certain disabilities.

The Interplak[®] by Conair[®] Cordless Portable Water Flossing System enhances oral hygiene and significantly contributes to healthy gums when used as part of a responsible oral health management program. It is particularly effective at cleaning around orthodontic appliances, fixed and removable bridge work, implants and splints, and it provides valuable gum massage with every use.

Remember: Oral irrigation is a supplement to, not a substitute for, your oral hygiene routine. Ask your dental health care provider to advise you on how the water jet can most effectively meet your specific oral hygiene needs.

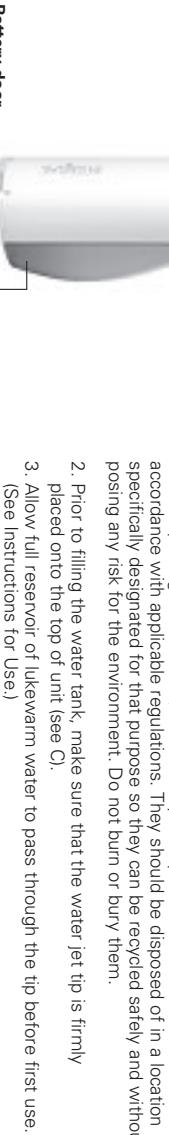Welcome to a New Level of At-Home Dental Care

The Interplak[®] by Conair[®] Cordless Portable Water Flossing System enhances oral hygiene and significantly contributes to healthy gums when used as part of a responsible oral health management program.

It is particularly effective at cleaning around orthodontic appliances, fixed and removable bridge work, implants and splints, and it provides valuable gum massage with every use.

Remember: Oral irrigation is a supplement to, not a substitute for, your oral hygiene routine. Ask your dental health care provider to advise you on how the water jet can most effectively meet your specific oral hygiene needs.

Initial Preparation:

Before using your Interplak[®] by Conair[®] Cordless Portable Water Flossing System for the first time:

1. Turn the battery door cover to the open position. Open door and insert 2 AA alkaline batteries (see A). Once installed, move battery door cover to the closed position.

Ensure that the batteries are installed correctly with regard to polarity (+ and -) (see B). Always purchase the correct size and grade of battery, most suitable for intended use. Replace all batteries of a set at the same time.

To avoid water getting into the battery compartment, the door must be firmly closed and **securely locked**. Grip the battery door knob and turn counterclockwise to tighten (see E).

NOTE: Use only alkaline batteries. Do not mix old and new batteries. Do not mix alkaline, standard (carbon-zinc), or rechargeable (nickel cadmium) batteries. Remove batteries from equipment that is not to be used for an extended period of time.

When replacing batteries, make sure you dispose of the old batteries properly. Clean the battery contacts and also those of the device prior to battery installation. Remove used batteries promptly.

When replacing the batteries, make sure you dispose of the old batteries in accordance with applicable regulations. They should be disposed of in a location specifically designated for that purpose so they can be recycled safely and without posing any risk for the environment. Do not burn or bury them.

2. Prior to filling the water tank, make sure that the water jet tip is firmly placed onto the top of unit (see C).
3. Allow full reservoir of lukewarm water to pass through the tip before first use. (See Instructions for Use.)

Instructions for Use:

1. Insert water jet tip onto the top of unit until you hear a "click" noise (see C). The water jet tip should not pull out easily. When changing tip, press the attachment release button and remove.

2. There is only one way to fill the water tank. Remove the water tank by sliding it down the unit (see D).





NOTE: You can also use your favorite dental rinse or diluted mouthwash. Thoroughly rinse reservoir with clean water after each use with mouthwash or dental rinse, and whenever unit has not been used for an extended period of time.

3. Begin with switch positioned to lowest setting (Lo) to allow gums to become accustomed to water jet cleaning. Before turning your unit on, place the jet tip in your mouth, keep your lips slightly closed, and **lean over the sink**. Let water flow from your mouth into the sink.

4. Adjust to higher setting (Hi) if desired, for more powerful water flow.

NOTE: If the pressure diminishes over time, replace batteries.
5. Turn handle to direct the water jet tip toward tooth surfaces and between teeth to rinse away plaque and trapped particles. The pulsing action of the water massages gums as it cleans the entire mouth.

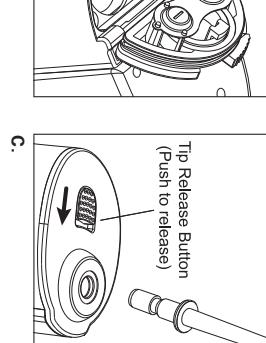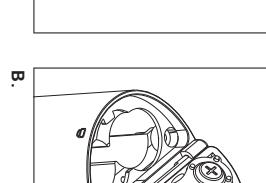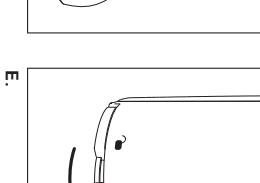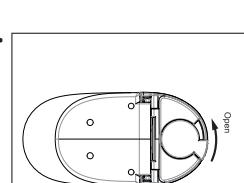6. If you travel, empty the water tank, switch off unit and take out the tip from unit. **NOTE:** To prevent damage to delicate gum areas, do not point water jet tip directly at gums. Water jet tip should be held $\frac{1}{4}$ in. away from your teeth and gums.

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