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CONSERVAR ESTAS INSTRUCCIONES

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OPERATING INSTRUCTIONS

DRY HEAT TREATMENT

This heating pad comes assembled (Figure 4).

This appliance has a polarized plug (one blade is wider than the other). As a safety feature, this plug will fit in a polarized outlet only one way. If the plug does not fit fully in the outlet, reverse the plug. If it still does not fit, contact a qualified electrician. Do not attempt to defeat this safety feature.

1. Plug the electrical cord into a 120V AC electrical outlet.
2. Slide the switch to ON. The power indicator light will illuminate to signal that the heating pad is on.
3. Place heating pad on top of the area to be treated.
4. When treatment is done, slide the switch to OFF and unplug. **NOTE:** This unit does not have an auto-off feature! Do not use or store heating pad in a folded position. Use and store lying flat.

- Consult your doctor before using in case of illness or medical condition, including but not limited to pacemaker use, pregnancy or cancer.
- If you experience pain in a muscle or joint for a prolonged period of time, discontinue use and consult your doctor. Persistent pain could be a symptom of a more serious condition.
- Never use any heating products on open wounds, discolored areas, or any area of the body that is swollen, burned, bruised, inflamed or where skin eruptions, ulcerated sores or varicose veins are present.
- This product should never be used by any individual suffering from any ailments that would limit or altogether inhibit the user's capacity to feel or have any sensation in the body part that this product is being used for.
- Use of this product should be pleasant and comfortable. Should pain or discomfort result, discontinue use and consult your physician.
- Ensure that your hands are dry when operating the controls or removing the plug.
- This product does not have an automatic off feature. Slide switch to OFF when done.

HOW TO ASSEMBLE

To assemble the heating pad, insert the heating pad completely into the cloth cover with the electric cord extending out from the cloth cover (Figure 4). If moist heat is desired, follow instructions in the Moist Heat Treatment section.


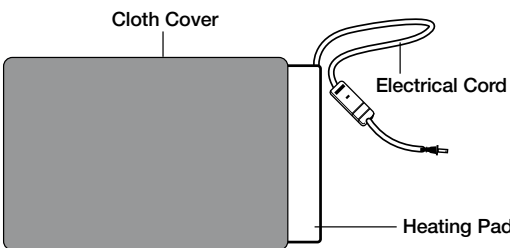
Figure 4

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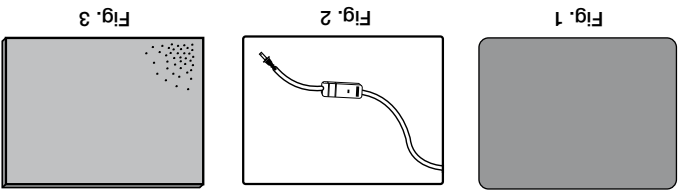SAVE THESE INSTRUCTIONS

HEATING PAD COMPONENTS

- Machine Washable Cloth Cover (Figure 1)
- Heating Pad (Figure 2)
- Absorbent Cloth (Figure 3)

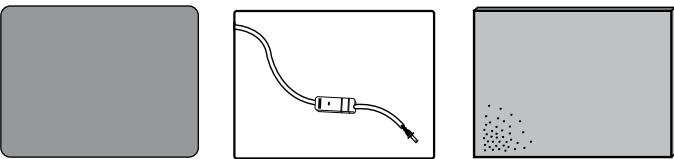
Figure 1

Figure 2

Figure 3

Cautions:

- The cloth cover must always be used with the heating pad. An absorbent cloth is provided for moist heat application.
- This heating pad should not be used in the following cases: If you have diabetes, tuberculosis, benign or malignant tumors, phlebitis and thrombosis, hemorrhages, inflammation of the skin, undetermined calf pain or joint pain, or if you have poor circulation. Consult a physician before using if you have any kind of fungal infection, swelling, fracture or persistent pain.
- **Consult your physician should you have any therapeutic questions regarding use of this heating pad.**
- Do not use on any unexplained pain, or swollen muscles, or following a serious muscular injury before consulting your physician.
- **Keep away from children.**
- Discontinue use of product if pain increases or swelling occurs after application.
- Do not use while bathing or showering. Use only in a dry area.

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TRATAMIENTO CON CALOR SECO

IMPORTANT SAFETY INSTRUCTIONS

DO NOT DESTROY

DANGER – To reduce the risk of burns, electric shock and fire, this product must be used in accordance with the following instructions:


1. READ ALL INSTRUCTIONS.
2. **ALWAYS UNPLUG THE UNIT FROM THE ELECTRICAL OUTLET IMMEDIATELY AFTER USING.**
3. DO NOT USE WHILE SLEEPING. AVOID FALLING ASLEEP WHILE USING THE PRODUCT.
4. DO NOT USE ON AN INFANT OR AN ANIMAL.
5. THIS PAD IS NOT TO BE USED ON OR BY ANY INDIVIDUAL WITH CERTAIN DISABILITIES, SLEEPING OR UNCONSCIOUS PERSON, OR A PERSON WITH POOR BLOOD CIRCULATION UNLESS CAREFULLY ATTENDED.
6. DO NOT USE ON AREAS OF INSENSITIVE SKIN.
7. BURNS CAN OCCUR REGARDLESS OF CONTROL SETTING; CHECK SKIN UNDER PAD FREQUENTLY.
8. NEVER USE PAD WITHOUT REMOVABLE COVER IN PLACE.
9. **DO NOT USE IN AREAS WHERE OXYGEN IS BEING ADMINISTERED.**
10. **DO NOT USE PINS OR OTHER METALLIC MEANS TO FASTEN THIS PAD IN PLACE.**
11. **DO NOT SIT ON OR CRUSH PAD – AVOID SHARP FOLDS.**
12. **NEVER PULL THIS PAD BY THE SUPPLY CORD.**
13. **DO NOT USE THE CORD AS A HANDLE.**
14. CAREFULLY EXAMINE INNER COVER AND LINE CORD BEFORE EACH

USE. DISCARD THE PAD IF INNER COVERING OR LINE CORD SHOW ANY SIGN OF DETERIORATION.

15. DO NOT TAMPER WITH THIS PAD IN ANY WAY. THERE ARE NO USER SERVICEABLE PARTS. IF, FOR ANY REASON, THIS PAD DOES NOT FUNCTION SATISFACTORILY, RETURN IT TO AUTHORIZED SERVICE CENTER AT ADDRESS GIVEN.
16. DO NOT USE THIS PAD WITH LINIMENT, SALVE OR OINTMENT PREPARATIONS THAT CONTAIN HEAT-PRODUCING INGREDIENTS, AS SKIN BURNS COULD RESULT. USE WHEREVER HOT APPLICATIONS ARE DESIRABLE FOR PERSONAL COMFORT AND WHEREVER RECOMMENDED BY YOUR PHYSICIAN OR DENTIST FOR THE RELIEF OF PAIN.
17. **LOOP CORD LOOSELY WHEN STORING; TIGHT WRAPPING MAY DAMAGE CORD AND INTERNAL PARTS.**
18. KEEP THE PAD AWAY FROM HEATED SURFACES.
19. PLACE PAD ON TOP AND NOT UNDER THE PART OF THE BODY NEEDING HEAT.
20. DO NOT PUT THIS PAD IN THE MICROWAVE.
21. DO NOT LEAVE PAD UNATTENDED, ESPECIALLY IF CHILDREN ARE PRESENT.
22. DO NOT USE PADS UNDERNEATH BEDCOVERS AND OTHER TYPES OF BLANKETS. DO NOT PLACE COVERS, SUCH AS SHEETS AND BLANKETS, ABOVE PADS. ONLY THE CLOTH COVER PROVIDED WITH THE PAD SHOULD BE USED.
23. DO NOT USE PAD ON SENSITIVE SKIN OR WHEN TAKING MEDICATION THAT CAN CAUSE SENSITIVE SKIN.

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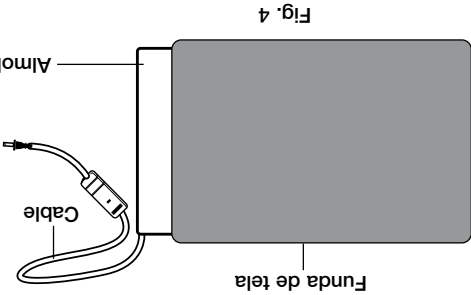
Indicador luminoso

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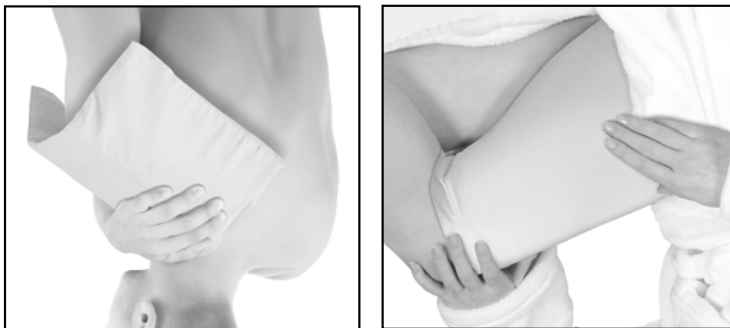NOTE: This heating pad does not shut off automatically. During use, check skin under heating pad frequently. Do not leave heating pad unattended. After use, dry the cloth cover, inner cloth and the heating pad before storing.

1. Ensure heating pad is unplugged.
2. Wet the absorbent cloth and squeeze out excess water.
3. Insert the damp cloth between the heating pad and the cloth cover.
4. Plug the electrical cord into a 120V AC electrical outlet.
5. Slide the switch to ON. The power indicator light will illuminate to signal that the heating pad is on.
6. Place heating pad on top of the area to be treated.
7. When treatment is done, press the ON/OFF button to turn unit off and unplug. For maximum benefit, use should not exceed 30 minutes at one time.
8. After use, dry the cloth cover, absorbent cloth and the heating pad before storing.

MOST HEAT TREATMENT

RESTRICTED HEATING PAD USES

For additional warnings, please see the Important Safety Instructions and Caution sections.

Not for use by children

Do not sleep while using heating pad

Do not use heating pad bunched up

Do not use under covers

Do not sit on heating pad

Do not sharply fold or crease heating pad

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LIMITED THREE YEAR WARRANTY (U.S. AND CANADA ONLY)

Conair will repair or replace (at our option) your unit free of charge for 36 months from the date of purchase if the appliance is defective in workmanship or materials.

To obtain service under this warranty, return the defective product to the service center listed together with your purchase receipt. California residents need only provide proof of purchase and should call 1-800-3-CONAIR for shipping instructions. In the absence of a purchase receipt, the warranty period shall be 36 months from the date of manufacture.

ANY IMPLIED WARRANTIES, OBLIGATIONS, OR LIABILITIES INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, SHALL BE LIMITED TO THE 36-MONTH DURATION OF THIS WRITTEN, LIMITED WARRANTY. Some states do not allow limitations on how long an implied warranty lasts, so the above limitations may not apply to you.

IN NO EVENT SHALL CONAIR BE LIABLE FOR ANY SPECIAL, INCIDENTAL, OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS OR IMPLIED, WHATSOEVER. Some states do not allow the exclusion or limitation of special, incidental, or consequential damages, so the above limitation may not apply to you. This warranty gives you specific legal rights, and you may also have other rights, which vary from state to state.

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Centro de servicio
Conair LLC
Service Department
7475 North Glen Harbor Blvd.
Glendale, AZ 85307

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21PA079983

IB-17383

Service Center
Conair LLC
Service Department
7475 North Glen Harbor Blvd.
Glendale, AZ 85307
Visit us on the web at www.conair.com
©2021 Conair LLC
Glendale, AZ 85307
stock image
21 PA079983

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For your safety and continued enjoyment of this product, always read the instruction book carefully before using.

instruction booklet for Models HP01XG

CONAIRcomfort™ standard heating pad

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