
IMPORTANT SAFETY INSTRUCTIONS DO NOT DESTROY

DANGER – To reduce the risk of burns, this product must be used in accordance with the following instructions:

1. DO NOT HEAT THIS PRODUCT WITH ANY MEANS OTHER THAN A MICROWAVE OVEN.
2. THE PAD COULD BE PERMANENTLY DAMAGED IF HEATED FOR MORE THAN THE RECOMMENDED GUIDELINE TIMES. THIS COULD AFFECT THE WARRANTY.
3. IF OVERHEATED, YOUR PACK WILL BE TOO HOT TO TOUCH. ALLOW IT TO COOL BEFORE REMOVING FROM THE MICROWAVE. CHECK IF THE PACK HAS BEEN DAMAGED. IF THERE IS ANY SIGN OF DAMAGE, DO NOT USE THE PRODUCT. IT SHOULD BE DISCARDED.
4. DO NOT USE WHILE SLEEPING. AVOID FALLING ASLEEP WHILE USING THE PRODUCT. DO NOT ALLOW CHILDREN TO USE OR HEAT THIS PRODUCT UNLESS UNDER ADULT SUPERVISION.
5. THIS PACK IS NOT TO BE USED ON OR BY ANY INDIVIDUAL WITH CERTAIN DISABILITIES, SLEEPING OR UNCONSCIOUS PERSON, OR A PERSON WITH POOR BLOOD CIRCULATION UNLESS CAREFULLY ATTENDED.

HOT/COLD BEAD PACK

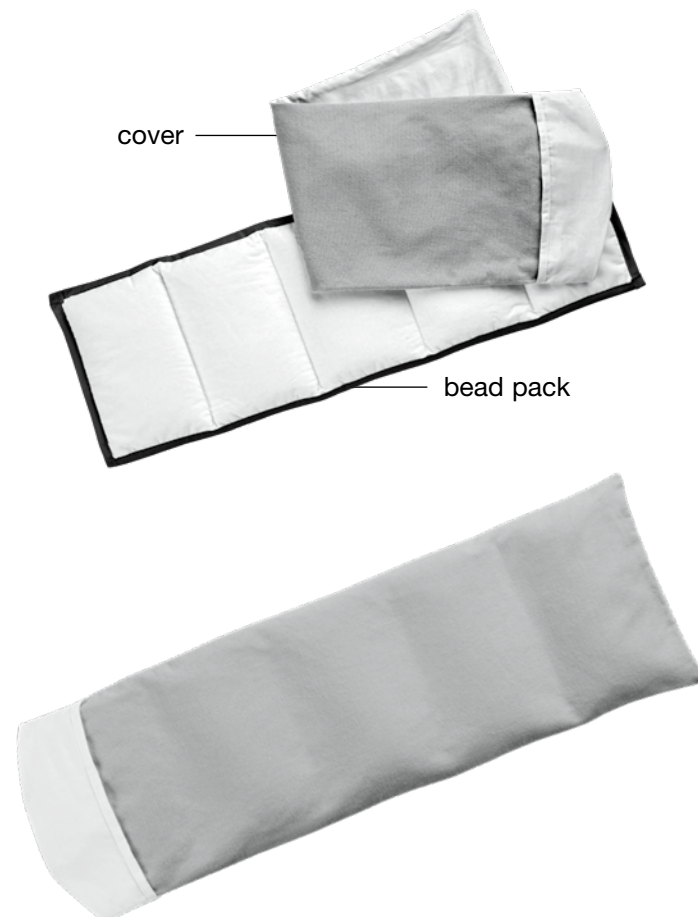6. DO NOT USE THIS PRODUCT WITH LINIMENT, SALVE OR OINTMENT PREPARATIONS THAT CONTAIN HEAT-PRODUCING INGREDIENTS. SKIN BURNS COULD RESULT. USE WHEREVER HOT APPLICATIONS ARE DESIRABLE FOR PERSONAL COMFORT AND WHEREVER RECOMMENDED BY YOUR PHYSICIAN OR DENTIST FOR THE RELIEF OF PAIN.
7. THIS PRODUCT SHOULD NOT BE USED BY DIABETICS, OR BY ANY PERSON WHO MAY BE INSENSITIVE TO THE APPLICATION OF DIRECT HEAT TO THE SKIN'S SURFACE, OR PERSONS WITH ANY SKIN DISORDER. IF IN DOUBT, CONSULT YOUR DOCTOR BEFORE USE.
8. THIS PRODUCT WILL PRODUCE MOISTURE. DO NOT USE IT RESTING AGAINST FURNITURE OR ANYTHING ELSE, AS THE PAD WILL MAKE THE SURFACE WET.
9. CAUTION: DO NOT EXCEED 30 MINUTES PER APPLICATION. RECOMMENDED TREATMENT TIME IS 10–30 MINUTES.
10. CHECK THE TREATMENT AREA EVERY 10 MINUTES TO HELP PREVENT BURNS OR FROSTBITE.
11. ALWAYS CONSULT YOUR PHYSICIAN BEFORE BEGINNING ANY HOT OR COLD THERAPY.
12. IF PAIN DOESN'T IMPROVE OR BECOMES MORE SEVERE OR IF NEW SYMPTOMS OCCUR, STOP USE AND CONSULT YOUR PHYSICIAN IMMEDIATELY.
13. DO NOT USE PINS OR OTHER METALLIC MEANS TO FASTEN THIS BEADED PACK IN PLACE.
14. CAREFULLY EXAMINE INNER COVER BEFORE EACH USE. DISCARD THE BEADED PACK IF INNER COVERING SHOWS ANY SIGN OF DETERIORATION.
15. DO NOT SIT ON OR CRUSH THE BEADED PACK.
16. NEVER USE THE BEADED PACK WITHOUT THE COVER IN PLACE.
17. READ AND FOLLOW ALL INSTRUCTIONS ON THE BOX OR THE INSTRUCTIONS PACKED WITH PRODUCT BEFORE USING.

HOT/COLD BEAD PACK

Pain relief for the whole body

Ceramic beads retain heat for 30 minutes and provide soothing therapy for the pain and stress of everyday life. Living a balanced life to maintain your health and well-being has become an important priority. The Hot/Cold Bead Pack is a natural way to ease the aches and pains of everyday life.

Unlike rice, oats, buckwheat or flaxseed, ceramic beads will not attract mold and mildew, will remain odor free and will not dry out or dissolve after multiple uses. The Hot/Cold Bead Pack delivers heat safely and effectively to provide soothing comfort.

Some suggested uses include:

- Upper back relief
- Shoulder relief
- Lower back relief
- Neck relief
- Knee relief
- Menstrual cramp relief

LIMITED THREE-YEAR WARRANTY (U.S. AND CANADA ONLY)

Conair will repair or replace (at our option) your unit free of charge for 36 months from the date of purchase if the appliance is defective in workmanship or materials.

To obtain service under this warranty, return the defective product to the service center listed, together with your purchase receipt. California residents need only provide proof of purchase and should call 1-800-3-CONAIR for shipping instructions. In the absence of a purchase receipt, the warranty period shall be 36 months from the date of manufacture.

ANY IMPLIED WARRANTIES, OBLIGATIONS, OR LIABILITIES, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, SHALL BE LIMITED IN DURATION TO THE 36-MONTH DURATION OF THIS WRITTEN, LIMITED WARRANTY. Some states do not allow limitations on how long an implied warranty lasts, so the above limitations may not apply to you.

IN NO EVENT SHALL CONAIR BE LIABLE FOR ANY SPECIAL, INCIDENTAL, OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS OR IMPLIED, WHATSOEVER. Some states do not allow the exclusion or limitation of special, incidental, or consequential damages, so the above limitation may not apply to you.

This warranty gives you specific legal rights, and you may also have other rights, which vary from state to state.

Service Center
Conair Corporation
7475 North Glen Harbor Blvd.
Glendale, AZ 85307

For information on any Conair product, call 1-800-3-CONAIR or visit us on the web at www.conair.com

Please register this product at www.conair.com/registration

©2020 Conair Corporation

20PA073071

SAVE THESE INSTRUCTIONS

Caution: Follow these instructions carefully. Misuse can cause severe burns. Overheating the ceramic bead pack in the microwave can damage or destroy the product.

Use Warm:

1. Ensure your microwave is clean and free of any food, debris, or liquids.
2. Insert the ceramic bead pack in the fabric cover and place in microwave.
3. In a 1,000+ watt microwave, heat the ceramic bead pack for 1 minute.
4. Remove from the microwave and massage the pad to distribute the ceramic beads evenly.
5. When heated, place the ceramic bead pack over the area in need.

Use Cool:

1. Place the ceramic bead pack in a sealed plastic bag; this will prevent any odors from being absorbed. (If stored in a sealed plastic bag, ceramic bead pack may be placed in the fabric cover.)
2. Store in refrigerator or freezer for a minimum of 1 hour.
3. When cooled, if not stored in the fabric cover, place the ceramic bead pack in the cover and place it over the area in need of treatment.

hot/cold bead pack

CONAIRcomfort™

Instruction Booklet for model HP46

For your safety and continued enjoyment of this product, always read the instruction book carefully before using.

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