## RECHARGEABLE POWER PLAQUE REMOVER

### Instruction Manual Model RTGX01 (All Versions)

For your safety and continued enjoyment of this product, always read the instruction book carefully before using.

Questions or Comments? 1-800-334-4031  
Visit us at our website: [www.conair.com](http://www.conair.com)

## IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

### READ ALL INSTRUCTIONS BEFORE USING

- DO NOT handle the plug with wet hands.
- DO NOT place or store the charger base where it can fall or be pulled into a tub, sink, or toilet.
- **DO NOT reach for a charger base that has fallen into water; unplug immediately. Do not reach into the water under any circumstances. Have the unit checked by an Authorized Service Center after it has been in water.**
- The unit is designed for use in a 120 VAC, 60 Hz outlet. Never force the plug into the outlet.
- DO NOT submerge the charger base in water or any other liquids.

### DANGER

To reduce the risk of electrocution:

- Close supervision is necessary when this product is used by, on, or near children or disabled persons.
- Never use while drowsy or while bathing or while sleeping.
- Keep the supply cord away from heated surfaces.
- Never operate the unit if it has a damaged plug, if it is not working properly, if it has been dropped or damaged, or if it has been immersed in water. Refer to an Authorized Service Center for repair.
- If the cord is damaged, discard the unit – the cord cannot be replaced.
- DO NOT excessively bend or pull the cord. Unplug the cord by grasping the plug – NOT THE CORD.
- Use this product only for its intended use as described in this manual. Only use attachments recommended by the manufacturer.
- This product is designed for teeth cleaning only – not for any other purpose.
- To prevent damaged bristles from breaking off while the unit is in use, replace your brush head if bristles are worn or bent.
- Never drop or insert any object into any opening or hose.
- Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.

**WARNING**  
To reduce the risk of burns, electrocution, fire, or injury:

## SAVE THESE INSTRUCTIONS

### MEDICAL WARNING

For your safety:

- Contact your dentist prior to use if you have had oral surgery in the previous twelve months.
- When you start any new oral hygiene program, your gums may initially bleed slightly. Contact your dentist if excessive bleeding occurs or if the bleeding continues after two weeks of use.
- Be sure that only the bristles touch your teeth and gums. Damage to teeth and gums may occur if any part of the unit touches them.
- This product is designed for oral cleaning only. DO NOT use the unit for any other purpose.
- As with any other new oral hygiene program, individuals with heart problems should check with their physician prior to using the RTGX01, and should always exercise particular care with any oral hygiene device.
- Avoid using a brush with crushed or bent bristles. Damaged bristles may break off while the unit is in use. In order to prevent this, replace your brush head every three to four months, or sooner if the bristles become worn or bent.
- The RTGX01 is a personal care unit and is intended for household use only.

### WARNING

Changes or modifications to this unit not expressly approved by Conair could void the user's limited warranty. Servicing should always be referred to an Authorized Service Center.

## Limited 2-Year Warranty (U.S. and Canada only)

Conair will repair or replace (at our option) your unit free of charge for 24 months from the date of purchase if the appliance is defective in workmanship or materials.

To obtain service under this warranty, return the defective product to the service center listed, together with your purchase receipt.

California residents need only provide proof of purchase and should call 1-800-334-4031 for shipping instructions. In the absence of a purchase receipt, the warranty period shall be 24 months from the date of manufacture.

**ANY IMPLIED WARRANTIES, OBLIGATIONS, OR LIABILITIES, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, SHALL BE LIMITED IN DURATION TO THE 24-MONTH DURATION OF THIS WRITTEN, LIMITED WARRANTY.** Some states do not allow limitations on how long an implied warranty lasts, so the above limitation may not apply to you.

**IN NO EVENT SHALL CONAIR BE LIABLE FOR ANY SPECIAL, INCIDENTAL, OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS OR IMPLIED, WHATSOEVER.** Some states do not allow the exclusion or limitation of special, incidental, or consequential damages, so the above limitation may not apply to you.

This warranty gives you specific legal rights, and you may also have other rights, which vary from state to state.

### Service Center

Conair LLC  
Service Department  
7475 N. Glen Harbor Blvd.  
Glendale, AZ 85307

Please register this product at: [www.conair.com/registration](http://www.conair.com/registration)



## Congratulations on your purchase of the Interplak® Rechargeable Power Plaque Remover

Get ready to experience a brand-new level of oral clean. State-of-the-art engineering has created a brush head with specially designed bristles that perform specific cleaning functions.

- Feathered subgingival bristles to clean under the gum line
- Contoured polishing tufts to smooth the tooth surface

With the push of one button, you get safe and effective three-way cleaning. Happy, healthy brushing, from Interplak!

### Advanced Brush Head Design
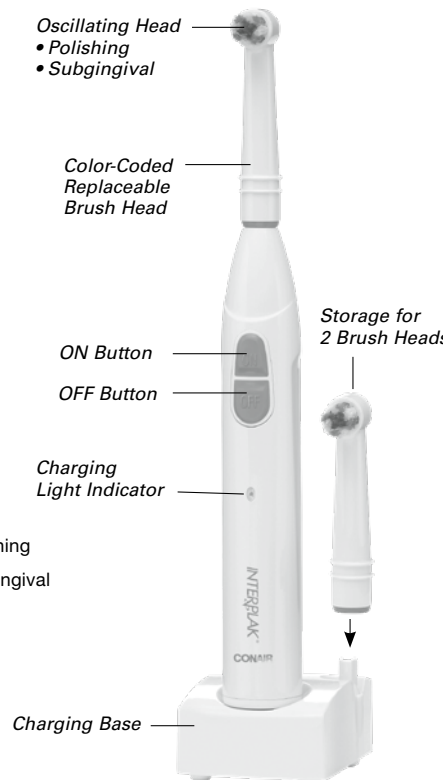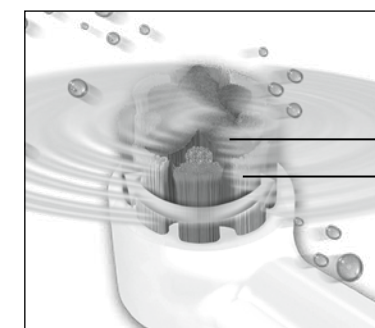
## Initial Preparation


Before using, charge the power handle for 24 hours. To charge, press the OFF button on the power handle and place it in charger base. The indicator light will go on.

NOTE: Be sure the electrical outlet is functioning. You can test it with another appliance. Also, check the electrical outlet reset button and/or wall switch.

About twice a year, completely run down and then recharge your power handle to maximize the life of the rechargeable battery.


## Instructions for Use: 5 Simple Steps

1. **Attach the brush head.**  
After the unit has been fully charged, grasp the brush head by the neck and place the power handle shaft inside the hollow end of the brush head; (A) press down and (B) turn clockwise until the brush head snaps into place. To remove, grasp the brush head by the neck (not the bristles), turn counterclockwise, and pull away from the handle. NOTE: It is recommended that you replace your brush head every three months, or sooner if the bristles become worn or bent.
2. **Moisten the bristles.**  
Hold the brush head under running water for a few seconds.
3. **Apply toothpaste sparingly and put brush head into mouth.**
4. **Turn the unit on.** Press the ON button.
5. **Position brush head properly.**  
Place bristles perpendicular to flat tooth surface and guide the unit slowly along the outside, then the inside, and finally the biting surfaces of the teeth. You do not have to press hard.
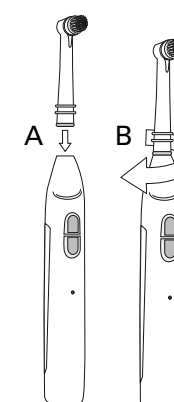
## Cleaning and Maintenance of Unit

1. **After each use**  
While toothbrush is on, hold the brush under running water for a few seconds.
2. **Turn the unit off**  
Press the OFF button.
3. **Towel-dry the power handle**  
Once dry, place the handle onto the plugged-in charger base to maintain full power.
4. **Brush head maintenance**  
Once a week, or anytime the power handle appears to be slowing down, turn power on and pour a small amount of mouthwash or water on the head.
5. This unit has no user-serviceable parts.
6. To store: Unplug the appliance from the outlet, store in appliance box or a clean, dry place. To avoid breakage, do not wrap the power cord around unit. Do not hang the unit by the cord.

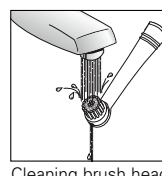Cleaning brush head

## Battery Safety Information

1. This product contains 1 nickel metal hydride (NiMH) AAA, 450 mAh 1.2V battery.
2. Since batteries are sensitive to high temperatures, when storing them, keep in a cool dry place.
3. Use only the power cord supplied with the appliance.
4. Never put batteries into a fire.
5. Switch off devices after use.
6. Keep batteries away from children, especially the batteries that are small enough to ingest.
7. If a battery is ingested, seek immediate medical help.
8. Never open, crush, puncture or disassemble batteries.
9. Dispose of batteries according to any State and Local regulations.
10. Do not dispose of rechargeable batteries as solid waste.

## Battery Disposal

1. As the batteries are rechargeable, you do not replace them. They are to be disposed of properly at the end of the toothbrush life. Please follow proper charging and disposal instructions provided.
2. To remove the battery for disposal or recycling, insert a coin into one set of slots on the battery end cap and twist counterclockwise to release the battery cap.
3. Firmly tap the bottom of the handle on a hard surface in order to slide out battery.
4. Now the battery is exposed, so you can separate the wire connectors. Carefully hold the wires on each side of the connector and gently cut the wires. The connector will separate and the battery can be removed from the inner housing. Bring the battery to the nearest collection center for disposal or call the toll-free RBRC® information line at 1-800-8-BATTERY. RBRC® is a registered trademark of the Call2Recycle organization. DO NOT throw the battery into the normal municipal waste system.

## FEDERAL COMMUNICATIONS COMMISSION NOTICE

**WARNING:** Changes or modifications to this unit not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment. This equipment has been tested and found to comply with the limits for an ultrasonic device, pursuant to Part 18 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. NOTE: This equipment has been tested and found to comply with the limits for an ultrasonic device, pursuant to Part 18 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

The device has been evaluated to meet general RF exposure requirement. To maintain compliance with FCC's RF exposure guidelines, this equipment should be installed and operated with a minimum distance of 0 cm between the radiator and your body.

## SUPPLIER'S DECLARATION OF CONFORMITY


Brand: INTERPLAK®  
Model No.: RTGX01  
Description: RECHARGEABLE TOOTHBRUSH  
Responsible Party: Conair LLC,  
1 Cummings Point Road, Stamford, CT 06902 1-800-334-4031  
Standards: FCC Part 18, ANSI C63.4-2014

This device complies with Part 18 of the FCC Rules. Operation is subject to the following two conditions:

- 1) This device may not cause harmful interference, and
- 2) This device must accept any interference received, including interference that may cause undesirable operation.

DATE OF ISSUE: XXXXXX

**CONTAINS NICKEL METAL HYDRIDE BATTERY. MUST BE COLLECTED, RECYCLED OR DISPOSED OF PROPERLY. DO NOT DISPOSE OF IN THE MUNICIPAL WASTE.**

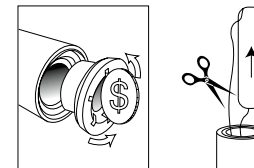For more information on battery replacement for your power toothbrush, call 1-800-334-4031 (8am to 5pm Mountain Standard Time, Mon-Fri).

Every 6 months, completely run down and recharge your power handle to maximize the life of the rechargeable battery.

To help protect the environment, this product is powered by a rechargeable nickel metal hydride battery.

Under various state and local laws, it is illegal to dispose of this battery in our municipal waste stream. The EPA-certified RBRC® Battery Recycling Seal on the nickel metal hydride (NiMH) battery indicates Conair is voluntarily participating in an industry program to collect and recycle these batteries at the end of their useful life, when taken out of service in the United States or Canada. The RBRC® program provides a convenient alternative to placing used NiMH batteries in the trash or the municipal waste stream, which may be illegal in your area. Please call 1-800-8-BATTERY for information on NiMH battery recycling and disposal bans/restrictions in your area. Conair's involvement in this program is part of our commitment to preserving our environment and conserving our natural resources or refer to the Call to Recycle link below for the nearest drop off location.

[www.call2recycle.org/locator/](http://www.call2recycle.org/locator/)

RBRC® is a registered trademark of the Call2Recycle organization.

## Power Brush Head Replacement

Dentists recommend replacing your toothbrush every 3 – 4 months, and this applies to replacing the brush heads of your Interplak® Rechargeable Power Plaque Remover. You can find 2 packs of Replacement Brush Heads, Model RBG3, at the same retail location where you purchased your power toothbrush. If you are having difficulty finding replacement brush heads, please call us at 1-800-334-4031 or visit us at [www.conair.com](http://www.conair.com).

## Product Features

Voltage: 120 VAC  
Power: 3W  
Hertz: 60 Hz

Battery: Nickel metal hydride rechargeable battery AAA type  
(explosion may occur if replaced with other type battery)

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